




Novaltis

*Seat/seat position transfer
Easy to handle and compact*

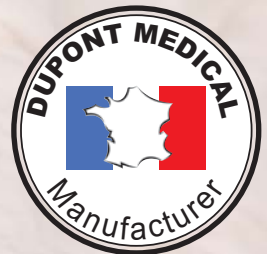

CE

Novaltis

Quick folding, discreet storage



No dismantling of the jack is necessary.
The jack remains protected by the structure of patient lift.

With its balanced and well-designed structure, it can be **vertically stored**, totally assembled and safely..

USEFUL AND PRACTICAL



POWER SUPPLY

- ① - Emergency stop
- ② - Electric emergency lowering
- ③ - Operating indicator



BATTERY CHARGING :

by connecting directly the device to a wall mains plug thanks to the integrated charger

or

at a distance by taking off the battery to connect it only in a wall mains plug.



Fast and effective transfer

Thanks to the new lumbar slings proposed on *Novaltis*, and to the pillow foam ample and deep holds-shins comfortable, the transfer will quickly be made and in complete safety.

This type of "semi-active" transfer particularly adapted to the patients having difficulty holding the standing position, presents advantages both for the transferred people and for the nursing :

- Stimulation of the residual autonomy.
- Beneficial lifting psychologically and physically, in particular for the transit.
- Change and toilet facilitated.
- Fast access and compactness of the product.



Light and easy to handle


Thanks to its minimum size and its strong and light structure, *Novaltis* is a manoeuvrable and secure patient lift.

Double-roller wheels,
75-mm Ø.
High-level rolling and
soundness.

Optimisation of the space at the base foot thanks to a
mechanical system without play.



Thanks to the ergonomic push-handle, people of all sizes can find their optimum hand position.

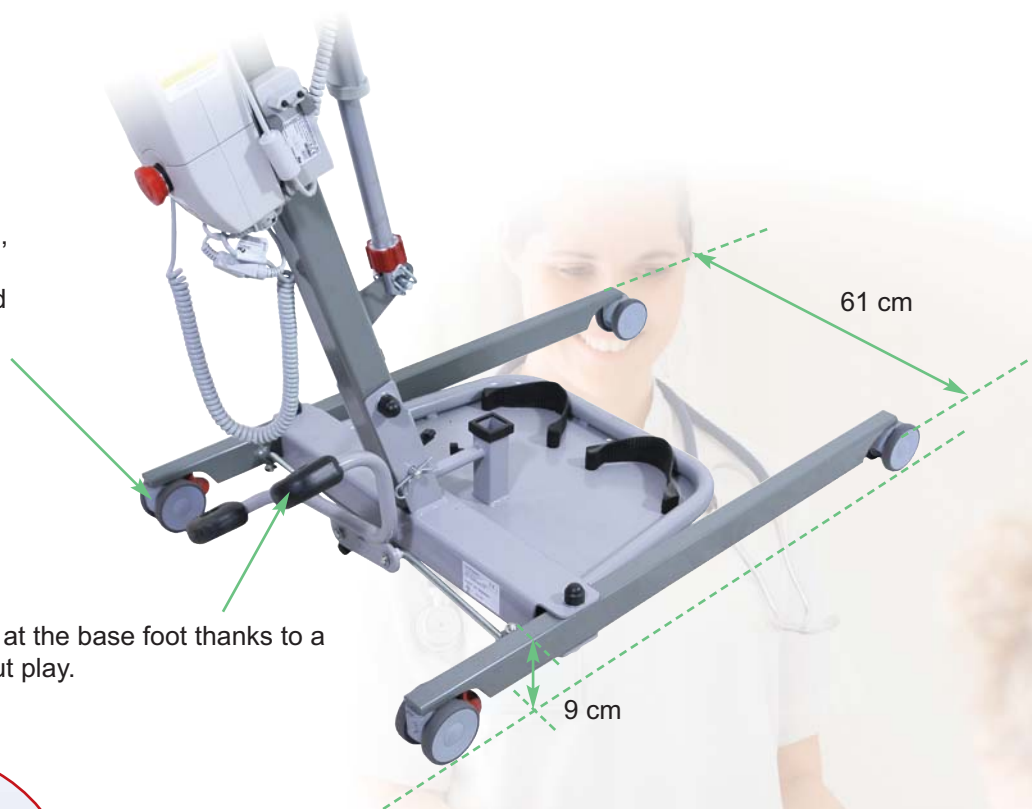

Two effective and ergonomic slings



Standard lumbar strap

Medium size - Réf. SR 370 1000

Large size - Réf. SR 370 2000



Lumbar strap with dorsal support

Medium size - Réf. SR 370 3000

Large size - Réf. SR 370 4000

SPECIFICATIONS

All the photographs are not contractual, and some components may change for manufacturing or improving purposes.

Reference of the Novaltis :	SR2910000
Lifting capacity :	175 kg
Cylinder :	8000 N
Overall width of the base :	61 cm
Overall length :	94 cm
Base height :	9 cm
Push handle height :	79 to 101 cm
Inner feet spacing :	54 cm to 81 cm
Patient lift weight :	39 kg



**On structure,
cylinder,
and control box**

Complies with NF - EN-standard / Cylinders IP X4 - 

DUPONT MEDICAL

Chaussée du Ban la Dame
Parc d'activités Eiffel Energie - ZAC du Ban la Dame
BP 19 - 54390 FROUARD (France)
Tel. +33 (0)3 83 495 495 - Fax. +33 (0)3 83 495 496
www.dupont-medical.com



SAMERY[®]
A DUPONT MEDICAL trademark